

# Darwin Chambers Company

http://www.darwinchambers.com

Phone: 314-771-3111

## KB074 ENVIRONMENTAL CHAMBER

Updated 2/05/05

### SPECS:

(AT 25°C / 50%RH)

Temperature Range	4° - 60°C
Temperature Control	±0.2°C
Temperature Uniformity	±1.0°C
Ambient Temperature Tolerance	17°C - 30°C
Control Sensitivity	0.1°C
Control Read Out	Actual and Set-Point Values
Temperature Sensor	NIST Traceable Platinum RTD
Controllers	Microprocessor Based PID
Calibration Capability	Yes, Actual and Set Values
Humidity Range	Ambient - 90% @ 25°C
Humidity Range with Optional Dehumidifica-	15-90% @ 25°C
Humidity Sensor	Vaisala NIST trace-
Control of Humidity at Sensor	±0.5%
Humidity Uniformity	±2.0%
Incubator Alarm	Audible and Visible
Alarm Type	Absolute Hi/Lo in
Remote Monitoring	Available
Electrical Requirements	115v, 60 Hz, 16 amps
Interior Cubic Feet	74
Required Clearance	1" All Sides

### CABINET CONSTRUCTION:

- 300 series stainless steel front grill, shroud and doors. Anodized quality aluminum exterior ends and back.
- Interior - INSF approved, white vinyl coated aluminum and 300 series stainless floor with coved corners.
- Cabinet mounted on a welded, galvanized frame rail, painted for additional corrosion protection. Frame rail fitted with 4" diameter casters, front casters include brakes.

### DOORS:

- 300 series stainless steel exterior with white aluminum interior to match cabinet liner. Doors extend full width of cabinet shell. Door gaskets of one piece construction. Glass doors available.
- Lifetime guaranteed adjustable torsion type closure system with aluminum block supports. The face of each door fitted with a 12" long recessed handle with sheetmetal interlock. Handles are foamed-in-place to insure permanent attachment. Barrel lock standard for each door. Locks keyed alike.

### INSULATION:

- Entire cabinet structure and doors insulated using foamed-in-place CFC free polyurethane foam.

### SHELVING:

- Adjustable vinyl coated wire shelves (3 per door). Shelves capable of supporting up to 250 lbs.
- Full height shelf support pilasters, with double oblong holes on 1/2" centers. Pilasters made of the same material as cabinet interior. Four (4) chrome plated, stainless steel shelf clips to be included for each shelf.

### LIGHTING:

- Incandescent interior lighting. Light(s) activated by rocker switch mounted above door(s).

### REFRIGERATION SYSTEM:



- Refrigeration system factory balanced using CFC free 134A refrigerant with top mounted condensing unit.
- Factory pre-engineered capillary tube system with oversized condenser and evaporator to reduce product pull down time and to consistently maintain interior cabinet temperature (non-adjustable capillary tube system superior to problematic expansion valve). Evaporator - epoxy coated at the time of manufacture to eliminate the potential of corrosion.
- Unsurpassed temperature control through continuous refrigeration bypass system.

### HUMIDITY

- Humidity control via microprocessor controlled ultrasonic humidification. Sensor is Vaisala, NIST traceable, w/cert.

### ELECTRICAL CHARACTERISTICS:

- Unit completely pre-wired at factory, and ready for final connection to a 115/60/1 phase - 20 amp dedicated outlet. A cord and plug set to be included.
- A sole use circuit is required.

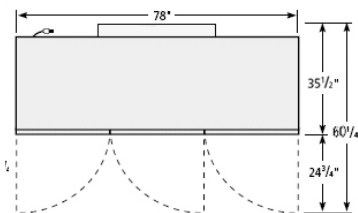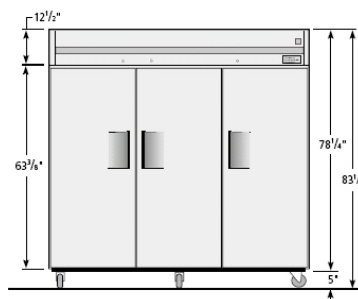

## ROUGH-IN DATA

Model #	Doors	Capacity	# of Shelves	Length	Depth	Height	Phase	Amps	Nema	Cord Length	Crated Weight
KB074	3	74 Cu. Ft.	9	78" 198cm	35.5" 90cm	83.3" 199cm	1	16	5-20	9'	550 lbs.